

ANTI-NEUROTROPHIN ANTIBODY SET 1

CATALOG NUMBER: NS112

LOT NUMBER:

QUANTITY: 500 µg each of AB1528SP, AB1513P, AB1517P, AB1519P and Negative Control

DESCRIPTION: The fate of many neurons in the developing and adult nervous system is controlled partly by neurotrophins. This family of proteins include Nerve Growth Factor (NGF), Brain Derived Neurotrophic Factor (BDNF), Neurotrophin 3 (NT3) and Neurotrophin 4/5 (NT4/5).

Neurotrophins are known to regulate the survival and differentiation of selective populations of neurons during embryonic development and maintenance of specific function in adulthood. The secretion of neurotrophins supports the selection of neurons during the thinning of neuronal population by programmed cell death (apoptosis). However, the roles of these neurotrophic factors in the brain and peripheral tissues are not very well understood. Therefore the availability of antibodies capable of neutralizing each neurotrophin will greatly aid further exploration of the action of these factors.

CHEMICON's Anti-Neurotrophin Antibody Set 1 consists of IgG fractions of sheep polyclonal antibodies for each of the four neurotrophins and a negative control and is conveniently packaged for use in biological inhibition analysis.

The Antibody Set 1 allows neutralization of neurotrophin activities for each of the neurotrophins. The levels of cross reactivity of each antibody for other members of the neurotrophin family is very low (Table 1).

Table 1. Neurotrophin antibody cross reactivities.

Antibody	ED ₅₀ (µg/mL)				
	NGF	BDNF	NT3	NT4	Ig
NGF	+	-	-	+/-	-
	(0.6)				
BDNF	-	+	-	-	-
		(10.5)			
NT3	-	-	+	-	-
			(5.5)		
NT4	-	-	-	+	-
				(3.6)	
Ig	-	-	-	-	-

Note: (-) represents no cross reactivity between antibody and neurotrophin (ED₅₀ > 200 µg/mL).
(+) represents antibody neutralization of neurotrophin.
(+/-) represents slight cross reactivity.

Sample protocol for using Chemicon's neurotrophin antibodies to inhibit biological activities

Preparation of single neurons


1. Dissect neural tissue to be studied (eg. dorsal root ganglia from E10 chick embryo) and place into calcium/magnesium free phosphate buffered isotonic saline (CMF-PBS).
2. Remove CMF-PBS.
3. Add fresh CMF-PBS containing 1% trypsin.
4. Incubate the tube at 37°C for 15-20 minutes and agitate at regular intervals.
5. Remove most of the liquid.
6. Wash tissue twice in Hanks containing 10% heat inactivated horse serum to inactivate residual trypsin.
7. Wash in CMF-PBS.
8. Remove liquid by pelleting tissue between washings via centrifugation at 2000xg for 1-2 min.
9. Dissociate tissue (trituration) by taking up the tissue in CMF-PBS with a siliconized pasture pipette heated to reduce tip diameter to a fine bore. Expel the contents while pressing the tip against the bottom of the tube.
10. Monitor trituration by examining a drop of the dissociated cells on a slide using an inverted phase contrast microscope.

Inhibition of neurotrophic activity with Chemicon's neutralizing antibodies

11. Preincubate for one hour at room temperature the neurotrophin at 10 ng/mL in BME with 1% fetal calf serum or 2 mg/mL of bovine serum albumin and anti-neurotrophin antibodies in increasing concentration (1-100 µg/mL is recommended).

Results can be plotted as shown below

The effect of anti-NT4 antibody on the neurite outgrowth of embryonic DRG promoted by NT4.

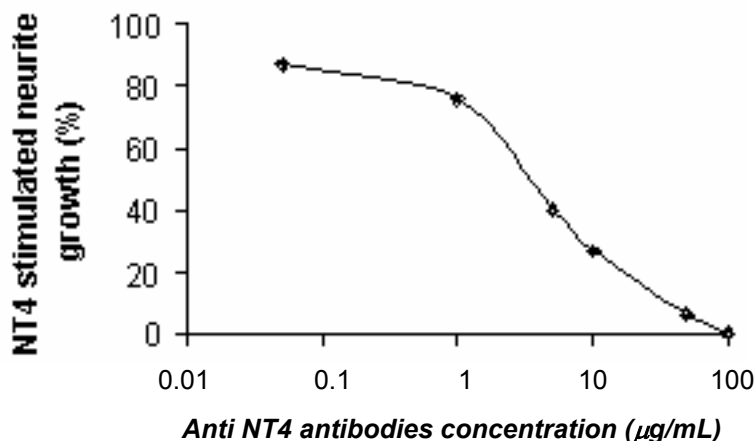

The ED50 of anti-NT4 neutralization is approximately 3.6 µg/mL

References:

Davis, A.M. (1989) Chapter 5: Neurotrophic factor bioassay using dissociated neurons. In Nerve Growth Factors, Edited by Rush, R.A. John Wiley and sons, New York.

Ebendal, T. (1989) Chapter 4: Use of collagen gels to bioassay nerve growth factor activity. In Nerve Growth Factors, Edited by Rush, R.A. John Wiley and sons, New York.

For research use only, not for use as a diagnostic.

Unless otherwise stated in our catalog or other company documentation accompanying the product(s), our products are intended for research use only and are not to be used for any other purpose, which includes but is not limited to, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses or any type of consumption or application to humans or animals.

**SHEEP ANTI-NERVE GROWTH FACTOR - BETA
(IgG Fraction)
POLYCLONAL ANTIBODY**

- CATALOG NO:** AB1528SP
- QUANTITY:** 500 µg
- SPECIFICITY:** Nerve growth factor-beta (NGF). By one site ELISA, less than 1% cross-reactivity against mouse NGF, recombinant human BDNF or neurotrophin 4. Approximately 2% cross-reactivity against neurotrophin 3.
- IMMUNOGEN:** Mouse nerve growth factor-beta purified from salivary glands. Antigen purity greater than 95% by PAGE.
- APPLICATIONS:** Immunohistochemistry: 1-10 µg/mL (protocol available).
Immunoblotting: 1-10 µg/mL
ELISA: 1-10 µg/mL
Inhibition of biological activity in vitro: 1-10 µg/mL
Optimal working dilutions must be determined by the end user.
- BIOLOGICAL ACTIVITY:** Neutralizes NGF, but not other neurotrophins.
- SPECIES REACTIVITIES:** Rat, chicken, mouse and human. Does not react with bovine NGF.
- FORMAT:** Protein G purified immunoglobulin.
- PRESENTATION:** Lyophilized from PBS, no preservatives. Reconstitute with 500 µL of sterile distilled water (1 mg/mL).
- STORAGE/HANDLING:** Maintain lyophilized material at -20 to -70°C for up to 12 months after date of receipt. After reconstitution maintain at -20 in undiluted aliquots for up to 6 months. Avoid repeated freeze/thaw cycles. Glycerol (ACS grade or better) can be added (1:1) for additional stability.
- REFERENCES:** *J. Neurosci. Meth.* **54**:95-102 (1994).

For research use only, not for use as a diagnostic.

Unless otherwise stated in our catalog or other company documentation accompanying the product(s), our products are intended for research use only and are not to be used for any other purpose, which includes but is not limited to, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses or any type of consumption or application to humans or animals.

**SHEEP ANTI-BRAIN DERIVED NEUROTROPHIC FACTOR
(IgG Fraction)
POLYCLONAL ANTIBODY**

CATALOG NO:	AB1513P
QUANTITY:	500 µg
SPECIFICITY:	Brain Derived Neurotrophic Factor (BDNF). By one site ELISA, less than 1% cross-reactivity against mouse NGF, recombinant human NT3 or neurotrophin 4.
IMMUNOGEN:	Recombinant Human Brain Derived Neurotrophic Factor (BDNF).
APPLICATIONS:	Immunohistochemistry: 1-10 µg/mL (protocol available). Immunoblotting: 1-10 µg/mL ELISA: 1-10 µg/mL Inhibition of biological activity <i>in vitro</i> : 1-10 µg/mL Optimal working dilutions must be determined by end user.
BIOLOGICAL ACTIVITY:	Neutralizes BDNF, but not other neurotrophins.
SPECIES REACTIVITIES:	Rat and human.
FORMAT:	Protein G purified immunoglobulin
PRESENTATION:	Lyophilized from PBS, no preservatives. Reconstitute with 500 µL of sterile distilled water (1 mg/mL).
STORAGE/HANDLING:	Maintain lyophilized material at -20 to -70°C for up to 12 months after date of receipt. After reconstitution maintain at -20 in undiluted aliquots for up to 6 months. Avoid repeated freeze/thaw cycles. Glycerol (ACS grade or better) can be added (1:1) for additional stability.

For research use only, not for use as a diagnostic.

**SHEEP ANTI-NEUROTROPHIN 3
(IgG Fraction)
POLYCLONAL ANTIBODY**

- CATALOG NO:** AB1517P
- QUANTITY:** 500 µg
- SPECIFICITY:** Neurotrophin 3 (NT3). By one site ELISA, less than 1% cross-reactivity against mouse NGF and recombinant human BDNF. Approximately 5% cross-reactivity against neurotrophin 4.
- IMMUNOGEN:** Recombinant Human Neurotrophin 3 (NT3).
- APPLICATIONS:** Immunohistochemistry: 1-10 µg/mL (protocol available).
Immunoblotting: 1-10 µg/mL
ELISA: 1-10 µg/mL
Inhibition of biological activity in vitro: 1-10 µg/mL
Optimal working dilutions must be determined by the end user.
- BIOLOGICAL ACTIVITY:** The antiserum will neutralize NT3, but not other neurotrophins.
- SPECIES REACTIVITIES:** Rat, chicken and human.
- FORMAT:** Protein G purified immunoglobulin
- PRESENTATION:** Lyophilized from PBS, no preservatives. Reconstitute with 500 µL of sterile distilled water (1 mg/mL).
- STORAGE/HANDLING:** Maintain lyophilized material at -20 to -70°C for up to 12 months after date of receipt. After reconstitution maintain at -20 in undiluted aliquots for up to 6 months. Avoid repeated freeze/thaw cycles. Glycerol (ACS grade or better) can be added (1:1) for additional stability.

For research use only, not for use as a diagnostic.

**SHEEP ANTI-NEUROTROPHIN 4
(IgG Fraction)
POLYCLONAL ANTIBODY**

CATALOG NO:	AB1519P
QUANTITY:	500 µg
SPECIFICITY:	Neurotrophin 4 (NT4). By one site ELISA, less than 1% cross-reactivity against mouse NGF and recombinant human BDNF. Approximately 5% cross-reactivity against neurotrophin 3.
IMMUNOGEN:	Recombinant Human Neurotrophin 4 (NT4).
APPLICATIONS:	Immunohistochemistry: 1-10 µg/mL (protocol available). Immunoblotting: 1-10 µg/mL ELISA: 1-10 µg/mL Inhibition of biological activity in vitro: 1-10 µg/mL Optimal working dilutions must be determined by the end user.
BIOLOGICAL ACTIVITY:	The antiserum will neutralize NT4, but not other neurotrophins.
SPECIES REACTIVITIES:	Human.
FORMAT:	Protein G purified immunoglobulin.
PRESENTATION:	Lyophilized from PBS, no preservatives. Reconstitute with 500 µL of sterile distilled water (1 mg/mL).
STORAGE/HANDLING:	Maintain lyophilized material at -20 to -70°C for up to 12 months after date of receipt. After reconstitution maintain at -20 in undiluted aliquots for up to 6 months. Avoid repeated freeze/thaw cycles. Glycerol (ACS grade or better) can be added (1:1) for additional stability.

For research use only, not for use as a diagnostic.

NEGATIVE CONTROL

CATALOG NO:	NS112-nc
LOT NUMBER:	xxxxxxx
QUANTITY:	500 µg
SPECIFICITY:	Sheep serum from an unimmunized animal.
APPLICATIONS:	For use as a negative control for the NS112 kit. Optimal working dilutions must be determined by the end user.
SPECIES REACTIVITIES:	No reactivity with rat. Other species have not yet been tested.
FORMAT:	Purified sheep immunoglobulin.
PRESENTATION:	Lyophilized, no preservatives. Reconstitute with 500 µL of sterile distilled water.
STORAGE/HANDLING:	Maintain lyophilized material at -20°C for up to 12 months after date of receipt. After reconstitution maintain at -20°C to -70°C in undiluted aliquots for up to 6 months. Avoid repeated freeze/thaw cycles. Glycerol (ACS grade or better) can be added (1:1) for additional stability.

For research use only, not for use as a diagnostic.

Unless otherwise stated in our catalog or other company documentation accompanying the product(s), our products are intended for research use only and are not to be used for any other purpose, which includes but is not limited to, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses or any type of consumption or application to humans or animals.