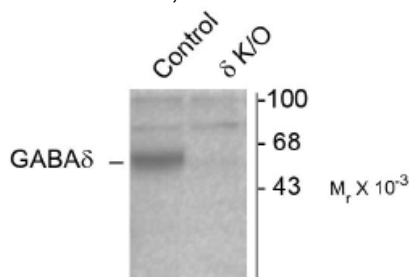


**RABBIT ANTI-GABA_A RECEPTOR, δ -subunit
AFFINITY PURIFIED
POLYCLONAL ANTIBODY****CATALOG NUMBER:** AB9752**LOT NUMBER:****QUANTITY:** 100 μ L**SPECIFICITY:** GABA_A Receptor δ -subunit. The antibody recognizes a protein of ~52 kDa in rat cerebellar lysates. The antibody did not react with a sample from a GABA_A Receptor δ -subunit knock-out animal.**IMMUNOGEN:** Fusion protein from the N-terminal of rat GABA_A Receptor δ -subunit.**APPLICATIONS:** Western blot: 1:1,000

Western Blot of rat cerebellar lysate. The immunolabeling is absent in lysate prepared from δ -knock-out animals.

Optimal working dilutions must be determined by the end user.

SPECIES REACTIVITIES: Rat, human, mouse and bovine. Other species have not been tested.**FORMAT:** Affinity purified immunoglobulin**PRESENTATION:** Liquid in 10 mM HEPES (pH 7.5), 150 mM NaCl, 100 μ g/mL BSA and 50% glycerol.**STORAGE/HANDLING:** Maintain at -20°C in undiluted for up to 6 months after date of receipt. Avoid repeated freeze/thaw cycles. Do not store in a self defrosting freezer.**RELATED REFERENCES:** Korpi ER, Mihalek RM, Sinkkonen ST, Hauer B, Hevers W, Homanics GE, Sieghart W, Luddens H (2002) Altered receptor subtypes in the forebrain of GABA_A receptor δ -subunit-deficient mice: recruitment of γ ₂-subunits. *Neurosci* 109:733-743.
McKernan RM, et al. (2000) Sedative but not anxiolytic properties of benzodiazepines are mediated by the GABA_A receptor α ₁-subtype. *Nature Neurosci* 3:587-592.
Mehta AK, Ticku MK (1998) Prevalence of the GABA_A receptor assemblies containing α ₁-subunit in the rat cerebellum and cerebral cortex as determined by immunoprecipitation: Lack of modulation by chronic ethanol administration. *Mol Brain Res* 67:194-199.

Ogris W, Pörtl A, Hauer B, Ernst M, Oberto A, Wulff P, Höger H, Wisden W, Sieghart W (2004) Affinity of various benzodiazepine site ligands in mice with a point mutation in the GABA_A receptor γ 2-subunit. *Biochem Pharmacol* 68:1621-1629.

Olsen RW, Tobin AJ (1990) Molecular biology of GABA_A receptors. *FASEB* 4:1469-1480.

Pörtl A, Hauer B, Fuchs K, Tretter V, Sieghart W (2003) Subunit composition and quantitative importance of GABA_A receptor subtypes in the cerebellum of mouse and rat. *J Neurochem* 87:1444-1455.

Wallner M, Hancher HJ, Olsen RW (2003) From The Cover: Ethanol enhances α 4 β 3 δ and α 6 β 3 δ γ -aminobutyric acid type A receptors at low concentrations known to affect humans. *Proc Natl Acad Sci (USA)* 100:15218-15223.

Whiting PJ, Bonnert TP, McKernan RM, Farrar S, Le Bourdellès B, Heavens RP, Smith DW, Hewson L, Rigby MR, Sirinathsinghji DJS, Thompson SA, Wafford KA (1999) Molecular and functional diversity of the expanding GABA_A receptor gene family. *Ann NY Acad Sci* 868:645-653.

Wohlfarth KM, Bianchi MT, Macdonald RL (2002) Enhanced Neurosteroid Potentiation of Ternary GABA_A Receptors Containing the δ -Subunit. *J Neurosci* 22:1541-1549.

Important Note: *During shipment, small volumes of product will occasionally become entrapped in the seal of the product vial. For products with volumes of 200 μ L or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the container's cap.*

For research use only; not for use as a diagnostic.

Unless otherwise stated in our catalog or other company documentation accompanying the product(s), our products are intended for research use only and are not to be used for any other purpose, which includes but is not limited to, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses or any type of consumption or application to humans or animals.

©2007: Millipore Corporation. All rights reserved. No part of these works may be reproduced in any form without permission in writing.