


**RABBIT ANTI-VAMP ASSOCIATED PROTEIN A
POLYCLONAL ANTIBODY**

CATALOG NUMBER: AB5741
LOT NUMBER:
QUANTITY: 200 µg
SPECIFICITY: VAMP associated protein A (VAPA).
IMMUNOGEN: Synthetic peptide from human VAMP associated protein A.
APPLICATIONS: Western blot: 1 µg/mL diluted in PBS containing 5% dry milk.

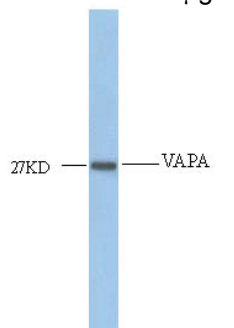

Hela whole cell lysate

ELISA, direct: 1:3,000
Optimal working dilutions must be determined by the end user.

SPECIES REACTIVITY: Human.
FORMAT: Purified immunoglobulin.
PRESENTATION: Lyophilized from PBS containing 2% sucrose. Reconstitute with 100 µL of sterile distilled water for a concentration of 2 mg/mL.
STORAGE/HANDLING: Maintain lyophilized material at 2-8°C for up to 6 months after date of receipt. After reconstitution maintain at -20°C in undiluted aliquots up to 6 months after date of receipt. Avoid repeated freeze/thaw cycles.

BACKGROUND REFERENCE: Weir, M.L. et al., *Biochem. J.* (1998) **333**:247-251.

Important Note: *During shipment, small volumes of product will occasionally become entrapped in the seal of the product vial. For products with volumes of 200 µL or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the container's cap.*

For research use only; not for use as a diagnostic.