


**RABBIT ANTI-VAMP ASSOCIATED PROTEIN B  
POLYCLONAL ANTIBODY**

**CATALOG NUMBER:** AB5665  
**LOT NUMBER:**  
**QUANTITY:** 200 µg  
**SPECIFICITY:** VAMP associated protein B (VAPB).  
**IMMUNOGEN:** Synthetic peptide from rat and mouse VAMP associated protein B.  
**APPLICATIONS:** Western blot: 1 µg/mL diluted in PBS containing 5% skim dry milk.

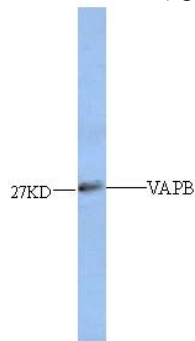

Mouse muscle extract

ELISA, direct: 1:3,000  
Optimal working dilutions must be determined by the end user.

**SPECIES REACTIVITY:** Rat and mouse.  
**FORMAT:** Purified immunoglobulin.  
**PRESENTATION:** Lyophilized from PBS containing 2% sucrose. Reconstitute with 100 µL of sterile distilled water for a concentration of 2 mg/mL.  
**STORAGE/HANDLING:** Maintain lyophilized material at 2-8°C for up to 6 months after date of receipt. After reconstitution maintain at -20°C in undiluted aliquots up to 6 months after date of receipt. Avoid repeated freeze/thaw cycles.  
**BACKGROUND REFERENCE:** Nishimura Y. et al., *Biochem. Biophys. Res. Commun.* (1999) **254**:21-26.

*For research use only; not for use as a diagnostic.*