


**CHICKEN ANTI-DNMT3a
AFFINITY PURIFIED
POLYCLONAL ANTIBODY**

CATALOG NUMBER: AB3431
LOT NUMBER:
QUANTITY: 100 µg
CONCENTRATION: 1 mg/mL
SPECIFICITY: DNMT3a (DNA methyltransferase 3 alpha isoform a)
IMMUNOGEN: Human DNMT3a, amino acids 1-300.
APPLICATIONS: Western blot: 1:500-1:2,000 (colormetric)

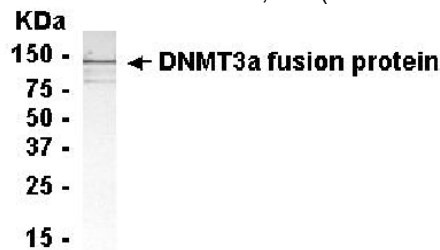

It is expected that the antibody will also work for immunohistochemistry and immunocytochemistry. Specific reactivity has not been confirmed. Optimal working dilutions must be determined by the end user.

SPECIES REACTIVITY: Human.
FORMAT: Affinity purified immunoglobulin.
PRESENTATION: Liquid in PBS.
STORAGE/HANDLING: Maintain at -20°C in undiluted aliquots up to 6 months after date of receipt. Avoid repeated freeze/thaw cycles.
BACKGROUND REFERENCE: Weisenberger, D.J., et al., *Gene* (2002) **298**:91-99.

For research use only; not for use as a diagnostic.

Important Note: During shipment, small volumes of product will occasionally become entrapped in the seal of the product vial. For products with volumes of 200 µL or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the container's cap.

**Now Available - Chicken IgY Negative Control Antibody (Cat. No. PP25),
Chicken IgY Blocking Solution (Cat. No. 2180) and Chicken IgY Precipitating Agarose (Cat. No. 2185)**