


CHEMICON

now part of Millipore

RABBIT ANTI-FBX2 AFFINITY PURIFIED POLYCLONAL ANTIBODY

CATALOG NUMBER: AB15184P
LOT NUMBER:
QUANTITY: 100 µg
CONCENTRATION: 1 mg/mL
SPECIFICITY: Fbx2. By Western blot the antibody recognizes a band at ~35 kDa on mouse brain lysate.
IMMUNOGEN: Synthetic peptide.
APPLICATIONS: Western blot: 0.5-1.0 µg/mL on mouse brain lysate.

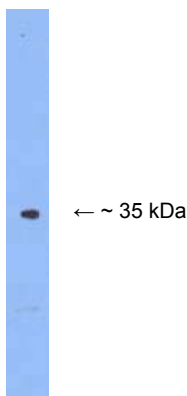

Optimal working dilutions must be determined by end user.

SPECIES REACTIVITIES: Mouse. Other species have not yet been tested.
FORMAT: Affinity purified immunoglobulin.
PRESENTATION: Liquid in PBS containing 0.1% sodium azide.
STORAGE/HANDLING: Maintain at 2-8°C in undiluted aliquots for up to 6 months after date of receipt.
IMPORTANT NOTE: During shipment, small volumes of product will occasionally become entrapped in the seal of the product vial. For products with volumes of 200 µL or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the container's cap.

For research use only; not for use as a diagnostic.

Unless otherwise stated in our catalog or other company documentation accompanying the product(s), our products are intended for research use only and are not to be used for any other purpose, which includes but is not limited to, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses or any type of consumption or application to humans or animals.

© 2007: CHEMICON International, Inc. - By CHEMICON International, Inc. All rights reserved. No part of these works may be reproduced in any form without permissions in writing.

USA & Canada • Phone: +1(800) 437-7500 • Fax: +1 (951) 676-9209 • Europe +44 (0) 23 8026 2233
Australia +61 3 9839 2000 • Germany +49-6192-207300 • ISO Registered worldwide
www.chemicon.com • custserv@chemicon.com • techserv@chemicon.com