

RABBIT ANTI-DCAMKL1 AFFINITY PURIFIED POLYCLONAL ANTIBODY

CATALOG NUMBER: AB15586

LOT NUMBER:


QUANTITY: 100 µg

CONCENTRATION: 1 mg/mL

SPECIFICITY: DCAMKL1 [Doublecortin-like and CAM kinase-like 1]. By Western blot the antibody recognizes a band at ~85 kDa on rat brain lysate.

IMMUNOGEN: Synthetic peptide.

APPLICATIONS: Western blot: 0.5-1.0 µg/mL on rat brain lysate.

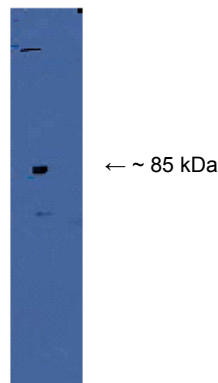

Optimal working dilutions must be determined by end user.

SPECIES REACTIVITIES: Rat. Other species have not yet been tested.

FORMAT: Affinity purified immunoglobulin.

PRESENTATION: Liquid in PBS with 0.1% sodium azide.

STORAGE/HANDLING: Maintain at 2-8°C in undiluted aliquots for up to 6 months after date of receipt.

IMPORTANT NOTE: During shipment, small volumes of product will occasionally become entrapped in the seal of the product vial. For products with volumes of 200 µL or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the container's cap.

For research use only; not for use as a diagnostic.

Unless otherwise stated in our catalog or other company documentation accompanying the product(s), our products are intended for research use only and are not to be used for any other purpose, which includes but is not limited to, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses or any type of consumption or application to humans or animals.